

**Sensible solution
for your patient**



AFIAS

Versatile Desktop Fluorescence Immuno-analyzer

AFIAS-6

- ▶ Automated POCT system with All-in-One cartridges
- ▶ Easy to use with proprietary C-tips (for fingertip blood)
- ▶ Up to 6 different parallel tests
- ▶ US FDA 510(K) cleared



Specifications

QC/Calibration

- 2 calibrators & 2 controls
- Internal & external controls

Interface

- 7" touch screen
- Built-in thermal printer
- USB port / Ethernet
- SD card slot / RS232
- LIS/HIS compatible

Memory

- 5,000 patient results
- 5,000 QC results
- 500 ID chips
- 100 User IDs

Analyzer

- 420 x 336 x 293 mm
- 15.1 kg
- 100 – 240V AC, 50-60Hz
- Internal temperature control

Automated All-in-one Fluorescence Immuno-analyzer

AFIAS-1

- ▶ Automated POCT system with All-in-One cartridges
- ▶ Easy to use with proprietary C-tips (for fingertip blood)
- ▶ Quick test results with reliability
- ▶ Small, compact, & light



Specifications

QC/Calibration

- 2 calibrators & 2 controls
- Internal & external controls

Interface

- 5" touch screen
- Built-in thermal printer
- USB port / Ethernet
- SD card slot / RS232
- LIS/HIS compatible

Memory

- 5,000 patient results
- 5,000 QC results
- 500 ID chips
- 100 User IDs

Analyzer

- 320 x 204 x 180 mm
- 3.9 kg
- 100 – 240V AC, 50-60Hz
- Internal temperature control

High-throughput Fluorescence Immunoassay with blood tube

AFIAS-10

- ▶ Full automation from blood tube
- ▶ 10 different tests in parallel
- ▶ Emergency test first during operation
- ▶ Bidirectional LIS/HIS



Specifications

QC/Calibration

- Internal & external controls

Interface

- 10.1" touch screen
- Built-in thermal printer
- USB port / Ethernet
- SD card slot / RS232
- LIS/HIS compatible

Memory

- 10,000 patients & QC results
- 100 ID chips
- 100 User IDs

Analyzer

- 426 x 443 x 395 mm
- 20 kgs
- 100-240V AC, 50-60 Hz
- Internal temperature control

Automated All-in-one Fluorescence Immuno-analyzer

AFIAS-3

- ▶ AFIAS system for 3 parallel tests
- ▶ Extra tip station – faster TAT
- ▶ Compact system for multiple testing
- ▶ User friendly UI



Specifications

QC/Calibration

- Internal & external controls

Interface

- 7" touch screen
- Built-in thermal printer
- USB port / Ethernet
- SD card slot / RS232
- LIS/HIS compatible

Memory

- 5,000 patients results
- 5,000 QC results
- 500 ID chips
- 100 User IDs

Analyzer

- 276 x 365 x 313 mm
- 9.3 kgs
- 100-240V AC, 50-60 Hz
- Internal temperature control

• All-in-One Cartridge

Ready-to-use reagent (single test format)

- ▶ Storage : 2-30 °C* for 20 months
- ▶ Sample type : Serum / Plasma / Whole blood
- ▶ Packing unit : 24 tests/Kit
- ▶ Test time : 3-22 min/each

* Some cartridges should be stored in 2-8 °C



• Can use whole blood, Finger prick or Venous

10 µL, 30 µL, whole blood

- No syringes
- No collection tube
- No pipettes
- No centrifuge

Flexibility of using C-tip™ as well as regular tube blood.




• 3 simple steps



Step 1: Collect blood




Step 2: Place the C-tip™

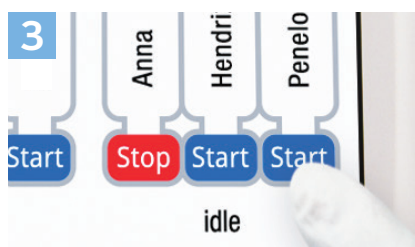

Step 3: Press START

• Study Results

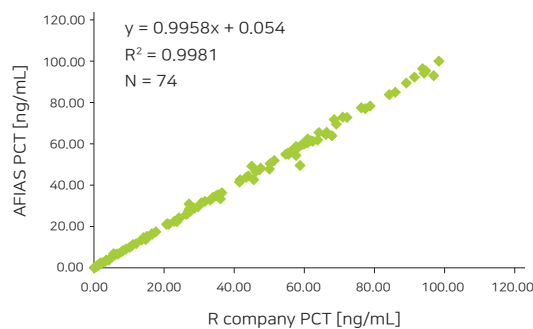

PCT Con. (ng/mL)	Intra-assay			Inter-assay		
	Avg.	SD	CV%	Avg.	SD	CV%
0.1	0.12	0.00	3.87	0.12	0.01	6.78
2.0	2.00	0.05	2.58	2.01	0.03	1.47
100	99.15	1.48	1.49	99.80	0.94	0.94

AFIAS Parameters

Category	Item	AFIAS-1 AFIAS-6	AFIAS-3	AFIAS-10	Sample type	Working range	Test time (min)
Cardiac	Tn-I Plus	●	●	●	WB/S/P	0.01-15.00 ng/mL	12
	CK-MB	●	●	●	WB/S/P	3-100 ng/mL	12
	CK-MB Neo	●	●	●	WB/S/P	3-100 ng/mL	12
	D-Dimer	●	●	●	WB/P	50-10,000 ng/mL	12
	D-Dimer Neo	●	●	●	WB/P	50-10,000 ng/mL	12
	NT-proBNP	●	●	●	WB/S/P	10-30,000 pg/mL	12
	Myoglobin	●	●	●	WB/S/P	2-500 ng/mL	12
	Myoglobin Neo	●	●	●	WB/S/P	2-500 ng/mL	12
	Cardiac Triple	●	●	●	WB/S/P	Tn-I: 0.01-15 ng/mL CK-MB: 3-100 ng/mL Myoglobin: 5-500 ng/mL	12
	ST2	●	●	●	WB/S/P	3.1-200 ng/mL	12
	BNP	●	●	●	WB/P	5-5,000 pg/mL	12
	Troponin T	●	●	●	WB/S/P	10-20,000 pg/mL	12
	hsCRP	●	●	●	WB/S/P	0.1-10.0 mg/L	3
Cancer	PSA	●			WB/S/P	0.5-100.0 ng/mL	15
	AFP	●	●	●	WB/S/P	5-350 ng/mL	15
	CEA	●	●	●	WB/S/P	1-500 ng/mL	15
Diabetes	HbA1c	●	●	●	WB	(NGSP) 4-15%	10
	HbA1c Neo	●	●	●	WB	(NGSP) 4-15%	10
	Microalbumin	●	●	●	Urine	2-300 mg/L	12
Hormone	TSH	●	●	●	WB/S/P	0.09-100 µIU/mL	12
	TSH Plus	●	●	●	WB/S/P	0.08-50 µIU/mL	12
	T3	●	●	●	S/P	0.77-7.7 nmol/L	10
	T4	●	●	●	S/P	10.23-300.00 nmol/L	10
	FSH	●	●	●	WB/S/P	1-100 mIU/mL	12
	Total βhCG	●	●	●	WB/S/P	5-50,000 mIU/mL	15
	β-HCG Plus	●	●	●	WB/S/P	2-1,500 mIU/mL	12
	LH	●	●	●	WB/S/P	1-100 mIU/mL	12
	PRL	●	●	●	WB/S/P	1-100 ng/mL	10
	Progesterone	●	●	●	S/P	4.45-127.2 nmol/L	15
	Testosterone	●	●	●	S/P	0.5-12.0 ng/mL	15
	Cortisol	●	●	●	WB/S/P	80-800 nmol/L	10
	AMH	●	●	●	S/P	0.02-15 ng/mL	12

WB = Whole blood, S = Serum, P = Plasma, NS = Nasopharyngeal swab

AFIAS Parameters

Category	Item	AFIAS-1 AFIAS-6	AFIAS-3	AFIAS-10	Sample type	Working range	Test time (min)
Infection	COVID-19 Ag	●	●	●	NS	Qualitative	12
	COVID-19 Ab	●	●	●	WB/S/P	Qualitative	10
	COVID-19/Flu Ag Combo	●	●	●	NS	Qualitative	20
	COVID-19 nAb	●	●	●	WB/S/P	Qualitative	20
	COVID-19 SP/NP IgG	●		●	WB/S/P	Qualitative	10
	CRP	●	●	●	WB/S/P	0.5-200 mg/L	3
	MxA/CRP	●	●	●	WB	MxA: 10-300 ng/mL CRP: 1-200 mg/L	12
	PCT	●	●	●	WB/S/P	0.1-100 ng/mL	12
	PCT X	●	●	●	WB/S/P	0.1-100 ng/mL	12
	PCT Plus	●	●	●	WB/S/P	0.02-50 ng/mL	12
	PCT Plus G	●	●	●	WB/S/P	0.02-50 ng/mL	12
	IL-6	●	●	●	WB/S/P	2-2,500 pg/mL	12
	Flu A+B	●	●	●	NS	Qualitative	10
	HBsAg	●	●	●	WB/S/P	Qualitative	12
	Anti-HBs	●	●	●	WB/S/P	Qualitative	15
	Anti-HCV	●	●	●	WB/S/P	Qualitative	12
	NORO	●	●	●	Feces	Qualitative	12
	ROTA	●	●	●	Feces	Qualitative	12
	Rota/Adeno	●	●	●	Feces	Qualitative	12
	IGRA-TB	●	●	●	P	Qualitative	15
Autoimmune	Anti-CCP Plus	●	●	●	WB/S/P	3.2-300 U/mL	12
	Total IgE	●	●	●	WB/S/P	1-1,000 IU/mL	12
TDM (Therapeutic Drug Monitoring)	Infliximab	●	●	●	WB/S/P	0.2-50 µg/mL	10
	Adalimumab	●	●	●	WB/S/P	0.2-50 µg/mL	10
	Golimumab	●	●	●	WB/S/P	0.2-30 µg/mL	10
	Trastuzumab	●	●	●	WB/S/P	3-150 µg/mL	10
	Bevacizumab	●	●	●	WB/S/P	5-300 µg/mL	8
	Total Anti-Infliximab	●	●	●	S/P	8-250 AU/mL	22
	Free Anti-Infliximab	●	●	●	WB/S	4-250 AU/mL	12
	Free Anti-Adalimumab	●	●	●	WB/S	3-200 AU/mL	12
	Free Anti-Golimumab	●	●	●	WB/S	2-40 ng/mL	10
	Free Anti-Trastuzumab	●	●	●	WB/S	8-500 AU/mL	10
Other test items	Ferritin	●	●	●	S/P	10-1,000 ng/mL	10
	Vitamin D	●	●	●	S/P	8-70 ng/mL	28
	Vitamin D Neo	●	●	●	WB/S/P	5-100 ng/mL	12

WB = Whole blood, S = Serum, P = Plasma, NS = Nasopharyngeal swab